

CRYOSKIN®

RÉVOLUTION



Cryo technology for non invasive,
non suction fat reduction,
body contouring and facelift

CRYOSLIM® | CRYOTONE® | CRYOLIFT® | CRYO CELLULITE | CRYO DOUBLE CHIN

Made in France

The potential of a natural technique

The business opportunity

Demand for aesthetic procedures is growing worldwide. Body reshaping and skin rejuvenation are among the most requested. However, going through the door of a plastic surgeon is still a barrier. This is why we believe that beauty salons and medspas are the best option to make aesthetic procedures available to everyone.

More than a cryolipolysis machine

Created in France in 2013, Cryoskin revolutionized the concept of cryolipolysis. Thanks to its handpiece using the Peltier technology, Cryoskin allows to destroy fat without aspiration and without damaging the skin while reshaping the body. In addition, Cryoskin offers body toning and facelift treatments, allowing for a faster return on investment.

Cryoskin Pact for Success



Face and
body



Attract
more customers



Rapid return of
investment

From evolution to Révolution


Cryoskin Révolution is a new concept. It is a modular platform that allows you to offer more treatments as your business grows.

The main unit allows to operate all Cryoskin treatments with the handpiece: fat reduction and body contouring, cellulite reduction, skin toning and lifting...

At any time you can add the static unit to offer hands-free fat removal treatments. It is not necessary to change the entire equipment, just connect the static unit to the main unit.

Cryoskin Révolution is manufactured in France, under CE standards, offering unmatched quality and technical performance.

Made in France

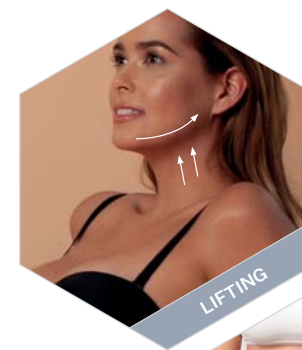

LIFTING



FAT REDUCTION




TIGHTENING



CELLULITE



DRAINAGE

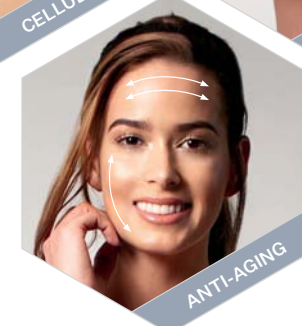

ANTI-AGING

Find out how Cryoskin Révolution can take you to a higher level...

The Cryoskin Method

CryoSlim®


CryoSlim® is a non-invasive procedure that eliminates fat and reshapes the figure at the same time.

CryoSlim® offers permanent results. Unlike other procedures that empty adipocytes and weaken the membrane, CryoSlim® completely destroys them, thus preventing them from reappearing.


FAT REMOVAL

A session lasts from 26 to 45 minutes depending on the area of the body and the program. There is an instant loss of centimeters at the end of the session, but the final result will be achieved 2 weeks later when the elimination of dead fat cells will be complete. It is recommended to do 5 to 10 sessions per body area depending on the amount of fat to be treated. Clients can lose up to 12 centimeters in 5 sessions.


A brief sequence of initial heat creates a dilation of the fat cells making them more receptive to the cold phase.

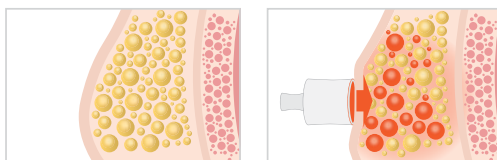

The second sequence begins with a rapid drop in temperature (Thermal Shock) leading to the apoptosis (fat cell destruction). Thermal Shock also stimulates lymphatic and blood exchanges.

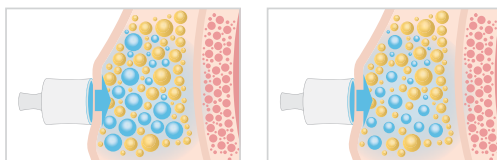

The final sequence uses high temperatures to reactivate the metabolism. The waste resulting from dead fat cells is disposed of through the lymphatic system and normal body processes within 2 weeks.

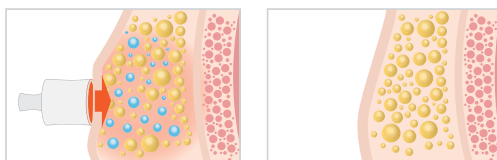


CLINICAL STUDY RESULTS

% of subjects showing improvement

overall satisfaction	94%
body shaping	87.5%
skin improvements	100%

Clinical study conducted by Dr Philippe Blanchemaison, MD, of the Société Française d'Accréditation Santé, Paris, France

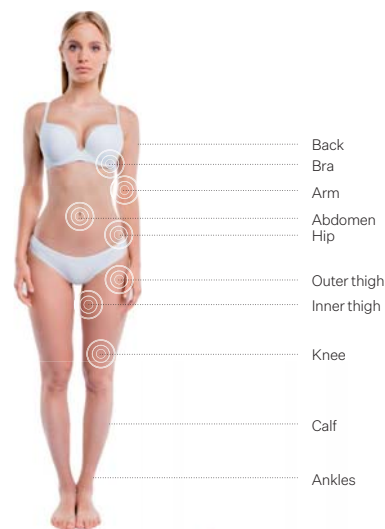

CryoSlim® abdomen



CryoSlim® hips



CRYOSKIN GUARANTEE FOR OPTIMAL RESULTS

- ✓ Thermal Shock
- ✓ Preset programs for each body area
- ✓ Exclusive protocols and massaging techniques



“

While trying Cryoskin myself, I lost an inch off each arm in the first session and that sold it to me straight away. Cryoskin impacted our business by a 30% increase in our revenue.

Ayse - Cryoskin owner, Enhance Clinic London

”


The Cryoskin Method

Cryo Cellulite

The Cryo Cellulite protocol combines slimming and toning techniques. Slimming eliminates the fat cells accumulated under the skin, responsible for the orange peel aspect. Toning accelerates lymphatic drainage, helping to eliminate fluids and toxins.


What is cellulite?

Cellulite appears when the adipocytes located under the hypodermis grow and group into clusters that exert pressure on the skin and lead to the formation of cellulite dimpling. At the same time, the hypodermis swells, compressing the blood and lymphatic vessels, preventing proper drainage.

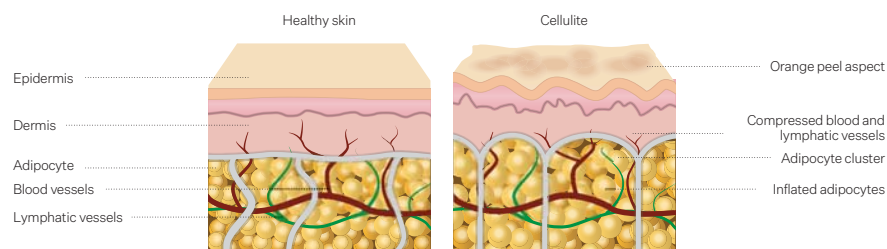


CRYO CELLULITE PROTOCOL

The slimming procedure reduces the size of adipocyte clusters, releasing the pressure on the skin. The skin appears smoother, showing an evident reduction of the orange peel aspect.

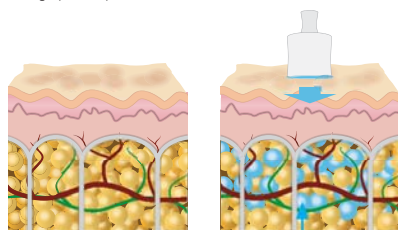

Apoptosis

The toning procedure accelerates lymphatic drainage in the treated area, draining out fluid and toxins. Blood supply increases toning up the skin and reducing edemas and swelling.

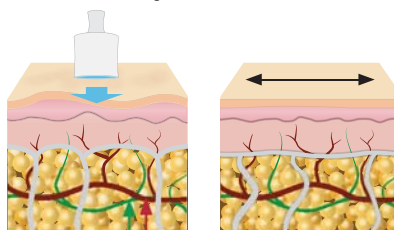

Increased blood and lymphatic flow

“

My regular clients are typically women in their 40s and 50s for whom cellulite is a major concern that is often associated with fluid retention.

What I like most about Cryoskin is the possibility of offering different treatments to meet the demands of each client.

Sandy - Cryoskin owner, Cryo Solihull

”

Cryo Cellulite



Before

After 1 session

Cryo Cellulite



Before

After 10 sessions

CLINICAL STUDY RESULTS

“This study demonstrated primarily a reshaping effect, in most volunteers. The reduction of the unappealing ‘orange peel’ aspect and the improvement of skin quality were significant.”

Clinical study conducted by Dr Philippe Blanchemaïson, MD, of the Société Française d'Accréditation Santé, Paris, France

CryoTone®

The CryoTone® uses low temperatures to tone the skin and improve local metabolism.

CryoTone® can be used for post-baby belly treatments, loose skin in the arms, cleavage etc. It is also possible to improve breast firmness. The cold has an instant tightening effect on the skin, but CryoTone® also acts in the long term by stimulating the production of collagen and elastin. It also performs a drainage of the treated area.



Before

After 8 session


The Cryoskin Method

CryoLift®

CryoLift® is a natural, non-invasive anti-aging solution. This procedure is particularly effective for tightening the skin of the face and neck, redefining the oval of the face. In addition, it provides a deep and lasting biological reaction by stimulating the natural production of collagen and elastin.

What is the skin aging process?

The aging process causes the skin's metabolism to slow down. Cells absorb fewer nutrients. Collagen and elastin fibers break down and the dermis loses elasticity, favoring the appearance of wrinkles. Additional factors such as sun damage and pollution can accelerate this process.

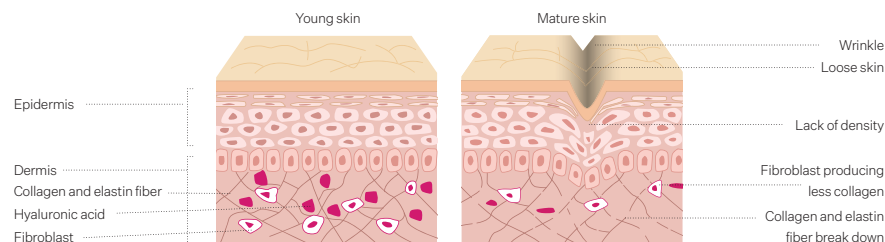


SKIN REJUVENATION

The CryoLift® technique creates a pumping effect on the blood and lymphatic systems, which increases oxygen and nutrient supply to the cells.

The fibroblasts benefit from higher nutrition and produce more collagen and elastin.

The dermis is re-densified and regains its elasticity. Wrinkles are naturally filled in. Skin surface is smoother and plumped up.

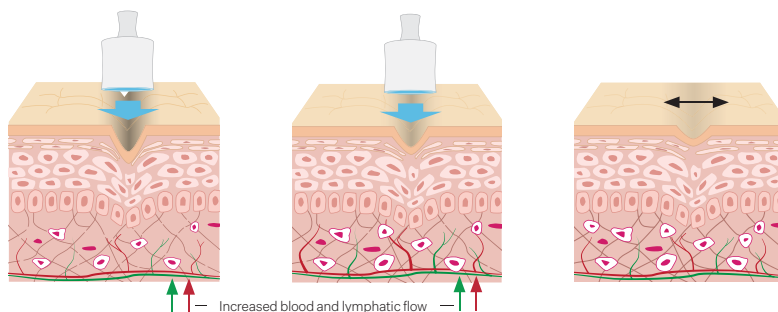

YOUNGER FOR LONGER

One session takes 20 minutes and the results are immediate. It is recommended to undergo one session a week for five weeks and then one session a month to maintain results. This plan has a cumulative effect that helps to maintain a higher production of collagen and elastin.


“

CryoLifts are a huge hit in our Medspa. The face is instantly lifted. It's a procedure I have once every two weeks at least. People are always willing to look after their skin, so the facial does provide consistency in revenue for us. With all the different procedures it provides, Cryoskin has been a great introduction to our business as it helped to increase our revenue.

Maria - Cryoskin Business Owner, LondonCryo

”

CryoLift®

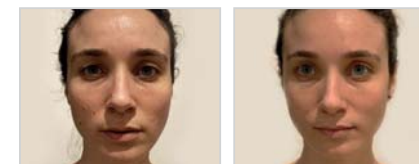

Before

After 1 session



Before

after 2 sessions

Cryo Double Chin


Cryo Double Chin is a specific protocol addressing excess of fat and loose skin in the neck area, redefining the oval of the face.

This method combines toning and slimming techniques, with 5 sessions of 15 minutes.



Before

After 5 sessions

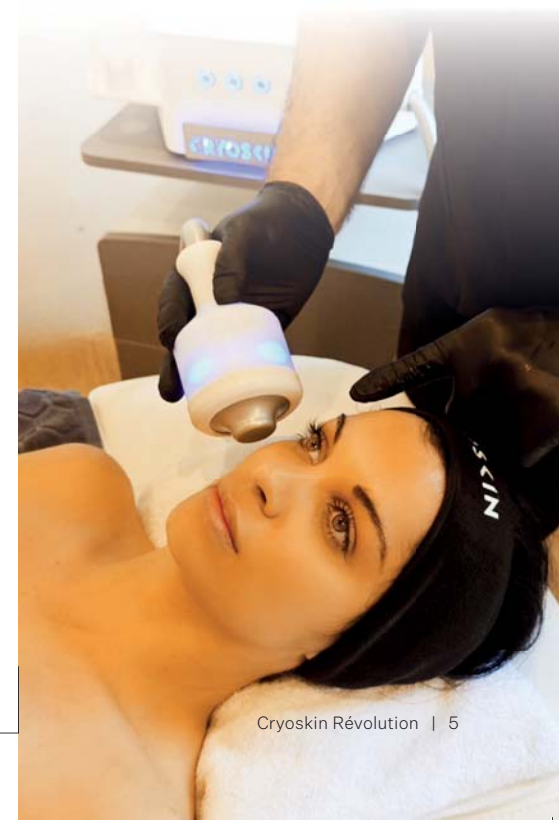

Main unit

Cryoskin Révolution

EFFICIENT

Cryoskin Révolution is the most powerful cryolipolysis machine on the market.

Thanks to a very efficient Peltier system, Cryoskin Révolution can deliver faster and stronger results.

SAFE

Made in France and complying with European CE quality standards, Cryoskin Révolution offers several security controls such as self-checks and automatic alerts in the event of a flow or temperature incident.

Made in France

COMFORTABLE

The programs are already preset for men and women, with different temperatures adapted to the different parts of the body. This allows to obtain the best results, being very comfortable for customers.

NEW ERGONOMIC HANDPIECE

- ✓ Light and compact
- ✓ Manufactured with aeronautical industrial hightech components
- ✓ Optional magnetic adaptor for facial treatments



NEW INTERACTIVE INTERFACE

- ✓ New software that allows you to choose the preset programs or select the Expert Mode to customize your own treatments
- ✓ Tutorial and videos integrated in the system
- ✓ Integrated customer database: query form, measurement record, etc.
- ✓ New WIFI connection that offers more features such as remote after-sales service and software updates



SPECIFICATIONS CRYOSKIN RÉVOLUTION MAIN UNIT

- Handpiece diameter: 65 mm / 2.56 inches
- Handpiece's cold transmission surface diameter: 55 mm / 2.16 inches
- Cord length: 2.5 m
- Unit overall dimensions: 44 x 40 x h21 cm including handpiece support
- Unit weight: 15.5 kgs
- Minimum temperature -10°C restrained to -8°C
- Maximum temperature +41°C restrained to +40°C
- Thermal Shock from +40°C to -8°C
- Unlimited memory to record data and new programs
- Dedicated cold generator
- Real time temperature sensor
- Electronic temperature management and temperature stabilizer
- Universal power supply: 110 – 230 V, 50/60 Hz
- Power consumption: maximum 300 VA
- LCD touch screen 10" powered by Android system
- Pre-set programs library
- Expert Mode to customize your programs
- Classification: Class 1 Electrical Appliance, Type B
- Made in France

Optional static unit

Offers hands-free treatments

The static unit allows to operate hands-free CryoSlim treatments on large body areas. You just have to place the 4 cryo heads on the body part to be treated and the machine will do everything by itself!

Depending on your type of clientele and your organization, you may decide to equip yourself with the static unit from the beginning or later. It is not necessary to change the whole machine, the units are completely compatible.

This flexibility is unique. Cryoskin Révolution adapts to your needs and your investment capacity.



COMBINED TREATMENTS FOR GREATER PROFITABILITY

- ✓ The 4 cryo heads can be used together or separately. This allows to treat two different areas
- ✓ They can also be used together with the handpiece, which allows to operate a manual CryoLift® together with a hands-free CryoSlim®



Easy to use pre-perforated extendable belts for hands-free procedures.

ENHANCED SAFETY

Allowing the operator to leave the client alone during hands-free procedures, without having to worry.

- ✓ New temperature stabilizer to avoid any variation during the hands-free sequences
- ✓ Stop hand button for the client, so they can stop the program in case of a problem when the operator is not in the room
- ✓ Restrained temperature range to avoid any risk of burn




Cryo heads vs suction cups

The suction cup technique was invented in the USA 30 years ago. It is the first generation of hands-free cryolipolysis devices. This technique often causes temporary or permanent damage to the skin, as it is strongly stretched within the cup. In fact, a period of 1 to 2 months is required between sessions to allow the skin to heal. Numerous clients who undergo suction cup procedures report experiencing pain.


With Cryoskin Révolution's hands-free mode, the 4 cryo heads remove fat without suction, stretching, burning or skin damage. The cold penetrates deep into the subcutaneous fat layers as the cryo heads are not moved during the entire treatment. Clients can undergo a session every 2 weeks (time needed to eliminate the dead fat cells). This higher frequency of visits allows for a greater turnover.

Suction cups

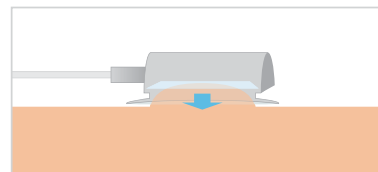

With suction cups, the area to be treated in one session is limited, since the cold can only reach the fat that is sucked inside the cup.

Cryoskin cryo heads

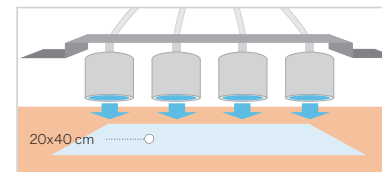


The 4 cryo heads radiate cold or heat under the skin, which makes it possible to cover areas of up to 20x40 cm / 8x16 inches, with a cold penetrating the subcutaneous fat layers up to 4 cm / 1.6 inch.

NEW TEMPERATURE STABILIZER

Other cold production devices work in on/off mode, which means that once the desired temperature is reached, the cooling stops. This produces temperature variations that can affect the result and, more importantly, expose the client to risk.

The Cryoskin Révolution temperature stabilizer allows to generate a stable temperature in each phase of the static treatment, which provides better results and greater safety.

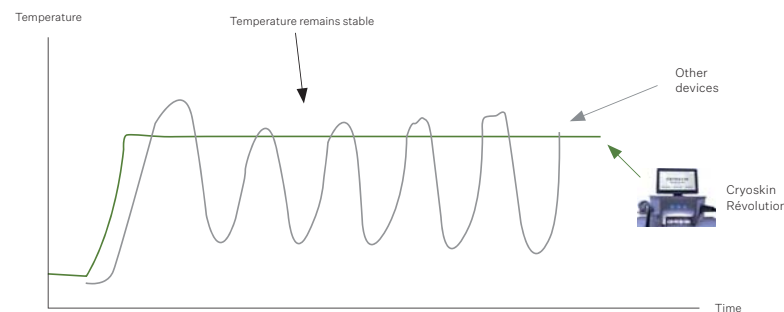

CRYOSKIN®

RÉVOLUTION

PRODUCT LINE GOING FROM DYNAMIC TO HANDS-FREE



CRYOSKIN DYNAMIC

Modular platform that allows you to offer more treatments as your business grows

- With the Dynamic handpiece, the beautician can perform all Cryoskin treatments: fat reduction, body contouring, cellulite reduction, skin toning and lifting treatments.
- At any time, you can upgrade your Cryoskin Dynamic and add the Cryoskin Static with a simple plug for dynamic, hands-free treatments.



CRYOSKIN STATIC

The hands-free unit for slimming treatments that saves time

- Place the 4 cryopads on the specific areas of the body that the customer wishes to treat and let the machine take care of everything!
- Each cryopad works individually, so you can use as many as you need and control the temperature of each one individually.
- At any time, you can upgrade your Cryoskin Static and add the Cryoskin Dynamic with a simple plug for dynamic, hands-free treatments.



CRYOSKIN STANDALONE

Standalone operation improves efficiency, allowing multiple treatments without constant supervision

- All Cryoskin treatments (dynamic and hands-free) in a single machine
- Integrated safety sensors prevent thermal shock and uneven temperature distribution.
- Automated, hands-free application reduces human error and improves treatment accuracy.
- Versatile technology supports cryolipolysis for body toning, skin tightening and improved circulation.



CRYOSKIN PARIS

@cryoskinparis

contact@cryoskinparis.com

+33 677626841

www.cryoskin-international.com

Cryoskin Révolution is manufactured in France, under CE standards, offering unmatched quality and technical performance.

Cryoskin Paris is part of the PRODESIGNPLUS group